User Manual



Important

Please read this User Manual, carefully before using the device or the service we provide.Failure to follow the instructions may seriously affect the performance of the product, or even cause damages.

Warnings

- The device complies with the relevant standards of electromagnetic compatibility.If you have medical concerns, please consult your doctor before using the product.
- Always use the cables and power adapter provided within the package.Otherwise it may not function properly or even cause damage to the product due to insufficient or excess power supply.
- Do not attempt to repair, disassemble, or modify the product without authorization.
- The device is not water resistant. Please avoid water and other liquid.

Content

1. About the Device 1

1.1. Inside the package 1

1.2. About the device 2

2. Product Installation and Configuration 3


3. Moving the device 5

4. Disclaimer 5


5. Warranty and Service 6

1. About the Device


1.1. Inside the package

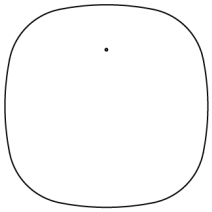


A

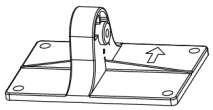


B

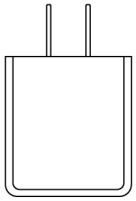


C

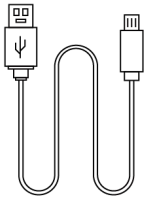


D

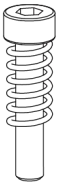

E

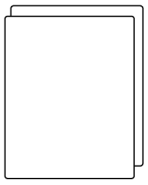


F

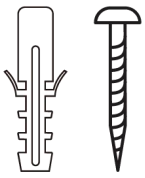

G

or

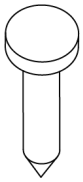

L

A: Sleep monitor

B:Mount

C: Power Adapter

D: Cord


E: Screw with Spring

F: Sticker


G:Screws and Wall plugs

L:Nail

1.2. About the device

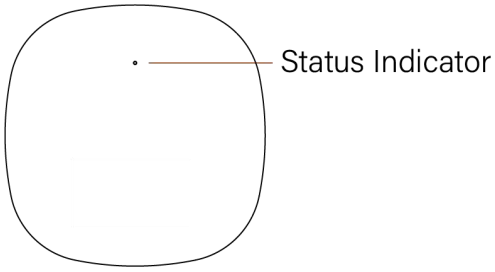

Device-FRONT

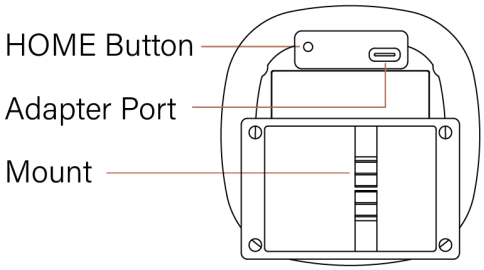

Device-BACK

Status Indicator	
Blue Light (flashing)	Device in setup mode.
Blue Light (constant)	Configuration is successful and the device is working properly.
Red Light (constant)	Device disconnected from the server or the Wi-Fi network.
Light Off	Once the device is successfully set up, the indicator will be blue for 30 seconds and then automatically turned off. The user can press the HOME button to wake up the indicator.

Product Parameters	
Product Model	See the label on the back of the device
Effective Detection Range	1.5m (4.92ft)
Communication Interface	2.4G Wi-Fi
Voltage Rating	DC 5V / 1A Max
Operating Temperature	-10-50°C (14-122°F)
Power Consumption	1.5W

2. Product Installation and Configuration

Tips

- Where and how to install the device are especially important. Improper installation may severely affect user experience. Be sure to follow the instructions for the best performance.
- The device supports only the 2.4GHz wireless networks. Please install it in a place where the Wi-Fi network signal is stable.

1) Where to install the device?

It is recommended to install the device in the bedroom, **right above the headboard and where the user to be monitored usually sleeps. Mount it on the wall, about 0.8-1.2 meter (2.6-3.9ft) away from the mattress**. (Fig1). The FRONT face of the device must aim directly at the user's chest, otherwise the result may not be accurate, The device is more accurate when there is only one user within the effective field of view.

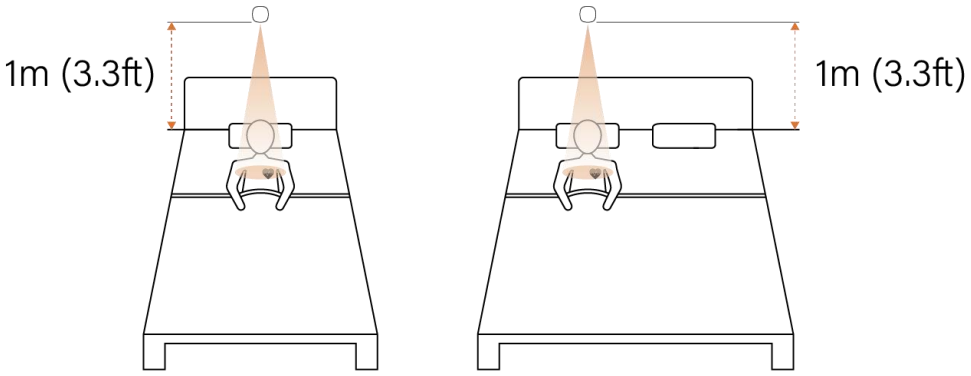


Fig.1

*It is recommended to install it about 1m (3.3ft) away from the mattress.

2) Secure the mount

Find a proper location to install the device. Please use the sticker and the nails to secure the mount on the wall. Please make sure the mount is horizontal and the arrow points up(Fig.2).

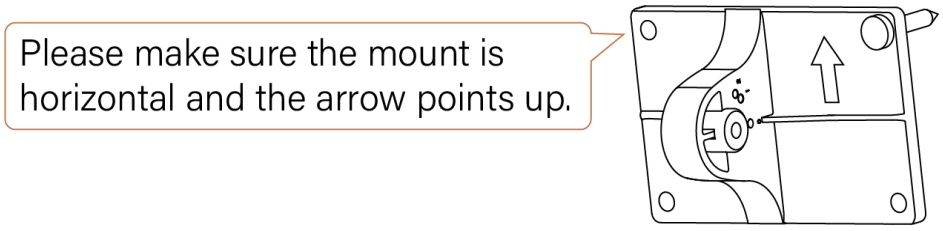


Fig.2

3) Device assembly

Fix the device to the mount using the screw with spring(Fig.3). Check whether the status indicator is on the top. If not, you may need to flip the mount.

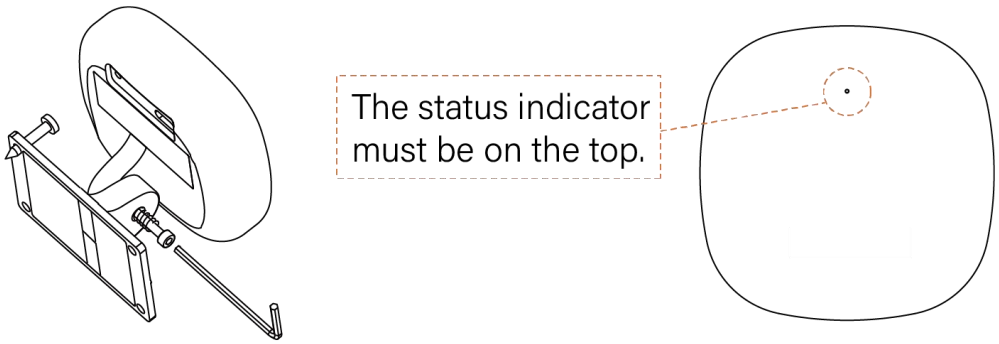


Fig.3

4) Adjust the device angle so that it aims directly at the chest of the user

For the best performance, please aim the device directly at the user's chest whilst he or she normally sleeps. To adjust the angle, please hold on to the device and push it to the right whilst tilting the device up or down to adjust the angle. If the device is installed around 08-1.2 m (2.6-3.9 ft), it is usually recommended to tilt the device at STOP1 Once the desired stop is reached the user will be able to feel that the locating pin locked in place.

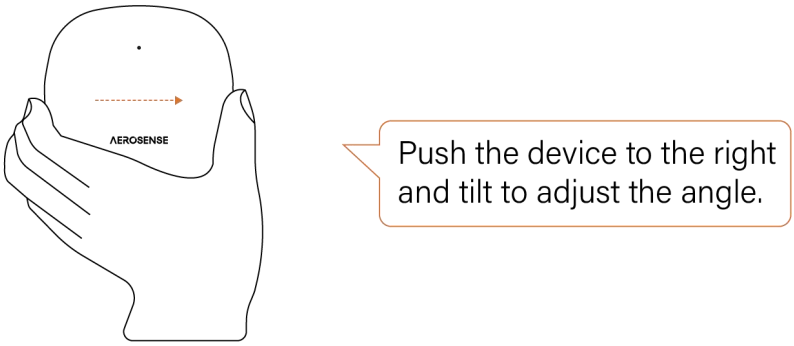


Fig.4

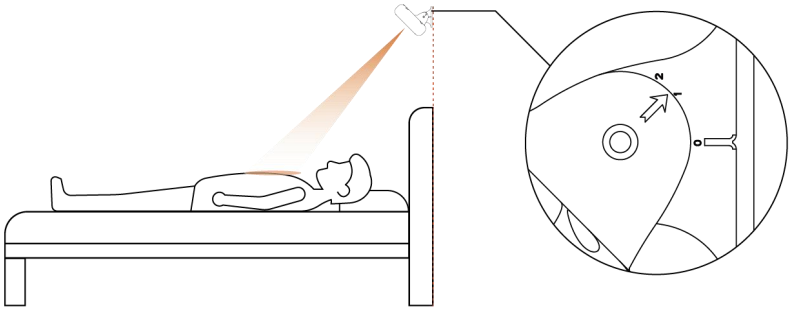


Fig.5

Adjust the angle so that the device is aiming directly at the chest of the user whilst the user is in his normal position whilst sleep.

* Please contact the sales or tech support for more information.

3. Moving the device

If the user needs to move the device to a different location, please follow the steps below:

1) Power off the device and remove the device and the mount.
2) Secure the mount and the device in the new location.
3) If the new location has the same Wi-Fi as the original one, the device will work automatically once powered back on; otherwise you will need to press and hold the Home button for 3-5 seconds until the blue light flashes and repeat the device configuration steps above.

4. Disclaimer

1) This product is not a medical device and does not guarantee 100% accuracy.
2) If the device is not installed properly, it will affect the accuracy and may result in

missed or false alerts.

3) The device may be less sensitive if the user is covered by a thick quilt, lies on his side or stomach.

4) The device will pick up all the vital signals within its field of view, Therefore pets and individuals other than the user who are present around the user may cause missed or false alerts.

5. Warranty and Service

When you purchase a new or certified refurbished device, we warrant the device against defects in materials and workmanship under ordinary consumer use for 1 year from the date of original purchase. The limited warranty applies only to hardware components of the devices that are not subjected to accidents, misuse, neglect, damages of any kind, alteration, repair, or commercial use. You can reach us through email or visit us at website.